




## Brachial Plexus Injury Rehab Support System

Neonatal brachial plexus injury is a complex orthopedic condition that often affects functioning alignment and joint development in the rest of the body, including the lower limbs and spine. The **Brachial Plexus Injury (BPI) Rehab Support System** is designed to extend the management of the affected arm and several associated postural mal-alignments into daily life for optimum training.

The TankTop provides a stable, Velcro®-sensitive surface area to attach strapping selected by the clinician to address such problems as postural deviations, scapular hypermobility, and shoulder subluxation. The Pelvic Anchor secures the TankTop against riding up, and increases the Velcro®-sensitive surface area for extending postural alignment strapping to the pelvis.

TheraTogs orthotic garment and strapping systems can also support and accelerate post-operative rehabilitation.

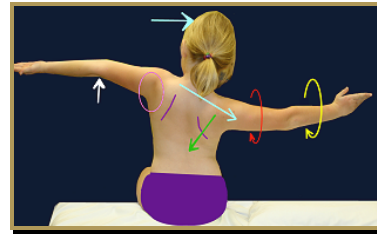

THE THERATOGS BPI SYSTEM CAN PROVIDE GENTLE, PROLONGED REPOSITIONING FORCES TO ADDRESS POSTURAL COMPENSATIONS FOR LIMITED ARM MOTION AND JOINT LAXITY.

If your client has these indications...	A BPI System can be used to...
Excessive movement in the torso when attempting to use the affected arm.	Reduce postural deviations to help to isolate the arm motion to the arm.
Flexible thoracic kyphosis	Maintain a more erect thoracic spine.
Pelvic obliquity with postural scoliosis due to muscle imbalances	Reduce lateral pelvic deviation. Check lower limb lengths to rule out the need for a shoe lift.
Shoulder protraction	Maintain improved shoulder alignment.
Scapular winging, elevation, and/or excessive abduction	Stabilize the scapula to help to isolate motion to the shoulder joint and promote soft tissue adaptation.
Inadequate scapulo-humeral disassociation	
Shoulder subluxation	Assist in holding the humerus in elevated position.

The BPI Rehab Support System includes components and instructions for both daytime functional retraining and protection for a lax shoulder joint during sleep.

## Model Sizes

Model No.	Description	Model No.	Description
#105	Infant BPI Rehab Support System: 15-25 lbs	#305	Pediatric BPI Rehab System: 45-65 lbs
#205	Pre-School BPI Rehab Support System: 25-45 lbs	#405	Petite Adult BPI Rehab System: 65-110 lbs

## How to Order

- ❑ In the USA, TheraTogs may be ordered through your local or regional orthotist or DME network, particularly if reimbursement is a requirement. See the TheraTogs website ([www.theratogs.com](http://www.theratogs.com)) for a list of Providers in your area.
- ❑ TheraTogs are also carried by leading distributors in the physical rehabilitation field all over the globe.

### In the USA:

FlagHouse • 800-793-7900  
[www.flaghouse.com](http://www.flaghouse.com)

Sammons Preston • 800-323-5547  
[www.sammonspreston.com](http://www.sammonspreston.com)

Abilitations / School Specialty • 800-850-8602  
[www.abilitations.com](http://www.abilitations.com)

TherAdapt • 800-261-4919  
[www.theradapt.com](http://www.theradapt.com)

### Outside the USA:

Argentina: Rehavita SA • Buenos Aires •  
[www.rehavita.com.ar](http://www.rehavita.com.ar)

Australia: Milestones Therapy • Sydney •  
[www.milestonestherapy.com.au](http://www.milestonestherapy.com.au)

Brazil: FisiVital • Marica • [www.fisiovital.com.br](http://www.fisiovital.com.br)

Canada: FlagHouse Canada • [www.flaghouse.ca](http://www.flaghouse.ca)

Canada: Sammons Preston Canada •  
[www.sammonspreston.ca](http://www.sammonspreston.ca)

Colombia: FisiMedica • Cali • [www.fisiomedica.com](http://www.fisiomedica.com)

Spain: Institut Sant Joan • Barcelona •  
[www.institutsantjoan.com](http://www.institutsantjoan.com)

India: Paediatric Physical Therapy Specialists Clinic •  
Narayanaguda, Hyderabad

Ireland / Scotland / UK: Rehab Logistics • Dublin •  
<http://www.rehab.ie/rehabenterprises/products.aspx>

Israel: Yaad Advanced Orthopedics, Ltd. • Tel Aviv •  
[www.yaadltd.com](http://www.yaadltd.com)

Italy: MOTUS SRL (formerly GOS Italia) • Verona •  
[www.gositalia.it](http://www.gositalia.it)

Poland: REH-RAD • Jelenia Góra • [reh-rad@wp.pl](mailto:reh-rad@wp.pl)

Venezuela: FisiMedica • Caracas • [www.fisiomedica.com](http://www.fisiomedica.com)

## About TheraTogs

TheraTogs provide a foam-lined, Latex<sup>®</sup>-free, breathable, comfortable, Velcro<sup>®</sup>-sensitive “external skin” for clinicians engaged in sensorimotor training to use as a field for attaching muscle-assist strapping for day-long “therapy you wear”.

TheraTogs are made of a patented, washable, composite fabric. The foam backing is made of an aqueous-based elastomeric urethane. TogRite™ strapping is Latex-free, elastomeric strapping with an inert, silicone-based grip surface.

TheraTogs™ are a Class I medical devices intended to be issued by, and applied under the supervision of, a licensed healthcare practitioner.

**TheraTogs are U.S. and foreign patent pending.**

**TheraTogs, Inc.** • 305 Society Drive, Suite C-3 • Telluride, CO 81435  
(888) 634-0495 • (970) 728-7028 fax • [www.theratogs.com](http://www.theratogs.com) • [info@theratogs.com](mailto:info@theratogs.com)